



The flexible and cost effective FM solution

With the Amigo 5 series, Oticon offers quality FM solutions in all price categories. The Amigo family is characterised by its reliability and ease of use, and the T5/R5 carry on this unique tradition.

Certainty The LED light on both the Amigo T5 and R5 provides instant visual confirmation that the transmitter and receiver are charged, and that an FM link has been established between the transmitter and the receiver.

Durability The internal components and external casing have been designed to withstand continuous use in demanding environments such as schools.

Flexibility The Amigo T5 and R5 are fully compatible with most FM receivers and transmitters, cochlear implants and hearing aids. By connecting headphones to the R5, children with **CAPD** (Central Auditory Processing Disorder), slight/mild hearing loss, and unilateral hearing loss can also benefit from using an FM system.

Convenience The operation and set-up of the Amigo T5 and R5 are extremely user friendly. The easy-to-use push buttons combined with a clear display make the system use and troubleshooting very easy.

Cost effective The T5 and R5 represent the most cost efficient solution in the Amigo family. So even more children can gain access to clear sound!


Amigo T5 transmitter

- External antenna built into the microphone cord
- One-touch mute
- Fail safe charging
- FM transmission range ~ 30 m



Amigo R5 receiver

- LED that confirms if the receiver is receiving an FM signal
- Headphones can be connected so that children with CAPD and mild hearing loss can benefit
- Easy connection to Cochlear Implants
- Multi-channel capabilities
- FM transmission range ~ 30 m
- SNR: 50 dB
- Fail safe charging

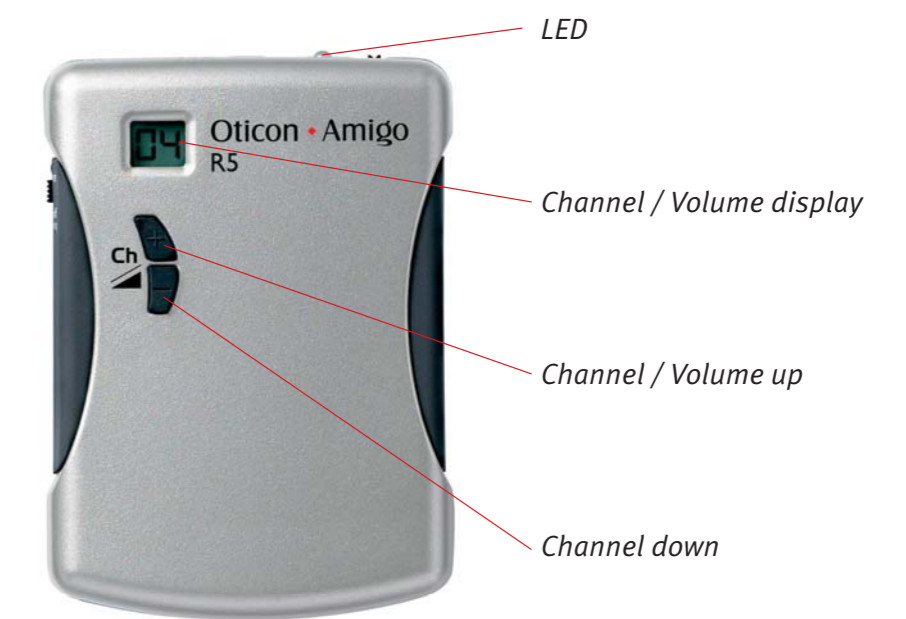

Amigo T5 has three different microphone options depending on your need:
- Directional lapel
- Omni directional lapel
- Boom



Amigo T5 and R5 come with a wall charger that can charge 2 units simultaneously.



The Amigo T5 can be connected to external audio sources (such as a TV, PC or MP3 player) with the aid of auxiliary cables.



Amigo R5 can be used with headphones for those diagnosed with CAPD.



The Amigo R5 offers an optional neck-loop for those using a Telecoil.



By connecting a DA1 cable to the Amigo R5, you can use it together with a hearing aid.

How Amigo can help

Tamper resistant volume and channel control

The Amigo R5 can be locked to prevent unwanted tampering with the volume and channel settings.

Fail-safe battery charging

The Amigo T5 and R5 use rechargeable AAA batteries, but it is safe to use standard alkaline batteries because of the intelligent battery detection scheme in the Amigo T5 and R5.

Quick & Easy Mute

The microphone can be muted when the FM signal is of no actual value - e.g. when the teacher is moving around in the classroom but not actually speaking.

Clear and simple user controls

The informative display, clearly marked jacks and easy to operate buttons make the Amigo system user friendly, even for inexperienced users.

LED on the Amigo

The built-in LED lights let you know that the signal is being received, the system is charging, the battery level is low or the transmitter is muted.

Children with CAPD

Children without hearing loss, but with CAPD, those with a slight/mild hearing loss or unilateral hearing loss often do not recognise subtle differences between sounds in words. By removing noise and transmitting the signal directly to their headphones, Amigo may improve their speech comprehension.



People First

We believe that it takes more than technology and audiology to create the best hearing instruments. That's why we put the individual needs and wishes of people with hearing loss first in our development of new hearing care solutions.

Oticon ♦ Amigo

20.01 / 00 01E 9H 016



FM made friendly
– and cost-effective