

**child
friendly
hearing
care**

It takes a truly dedicated approach to help children with hearing problems achieve their full potential. That's why we deliver the solutions and services that professionals and caregivers need to provide children the opportunities they deserve. This is what child-friendly hearing care is all about.



Amigo Arc

child
friendly
hearing
care

More of what
matters in FM

A choice that matters



What is Amigo Arc?

Amigo Arc is a universal neckloop receiver that acts as a bridge between an FM transmitter and the hearing aid. Using the neckloop it works with any instrument with a telecoil, such as BTE, RITE, ITE, cochlear implants or bone anchored instruments.

Why Amigo Arc?

Although FM has proven to be beneficial in classrooms, connecting these systems to hearing aids is not always easy. Amigo Arc, however, is simple and intuitive to use, offering consistent access to speech, day after day.

Higher bandwidth - more natural sound

For speech intelligibility, nothing beats Amigo Arc. Its extended, 7000 Hz bandwidth captures far more of the important HF speech cues than neckloops with a range of only 5500 Hz.

More options - higher flexibility

Amigo Arc offers superior sound quality not only for FM but also for music. It has a unique AUX IN option, where the bandwidth increases to 9000 Hz - plus a bass boost that promises far richer sound quality than any other neckloop currently available.

Another exclusive feature is the selection of colourful, exchangeable covers that offer children a chance to put their own personal touch on Amigo Arc.

The Amigo family at a glance



Amigo Arc



Amigo T21/T20: Advanced, robust FM transmitter and channel synchroniser in one. With Amigo T21 a built-in team teaching function is available.



Amigo T10: This elegant, multi-purpose FM transmitter offers three microphone modes: Omnidirectional, Focus and Super Focus.



Amigo T5: This sensible, easy-to-use transmitter fulfils all the basic FM requirements.



Amigo WRP: Optional accessory for programming and synchronising Amigo Arc.




Amigo eZync: This wall-mounted channel synchroniser automatically synchronises Amigo Arc and other Amigo receivers as students pass by.



Black Dark red Dark blue White Silver

Personalise Amigo Arc with colourful shells that are easy to click on and off.

A universal solution

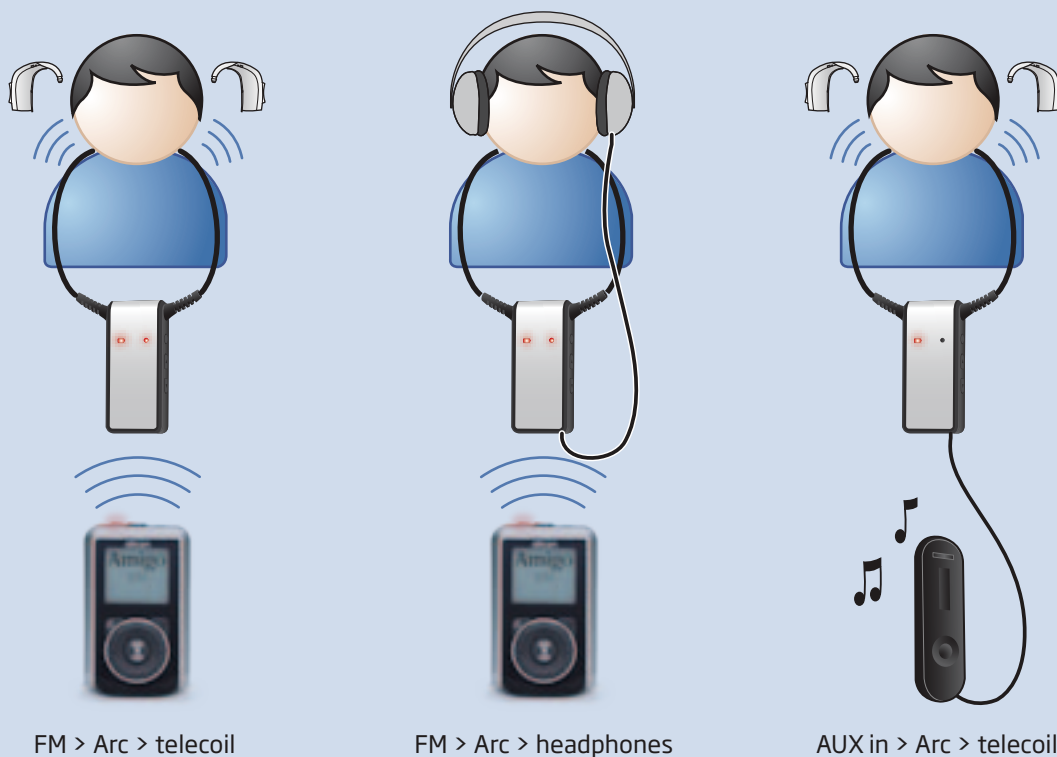



Making life easier at school

By acting as a bridge between the FM transmitter and the telecoil instrument, Amigo Arc offers children a significant advantage in class. The extra boost can also aid concentration in children with normal hearing, e.g. children with auditory processing disorder or attention difficulties.

A very useful feature in Amigo Arc is the visual status indicator, which provides instant confirmation that the FM system is working and the child is hearing clearly. This makes both caregiver and child feel more confident.

Making life at home more entertaining

Whenever children want to relax at home they can use Amigo Arc in one of two ways: they can connect it directly to a music player, or they can pick up the signal from an FM transmitter that's connected to their TV. This allows them to listen at the same volume as the rest of the family.

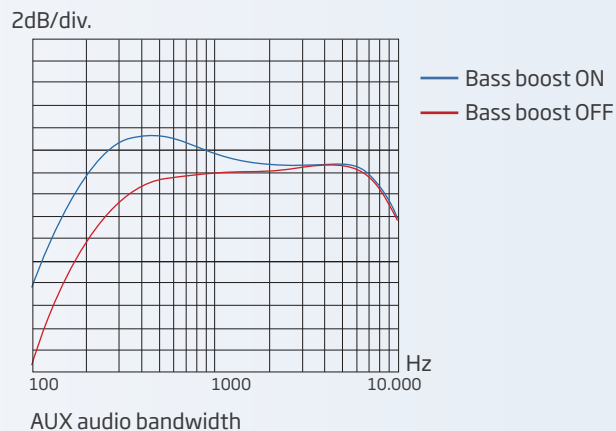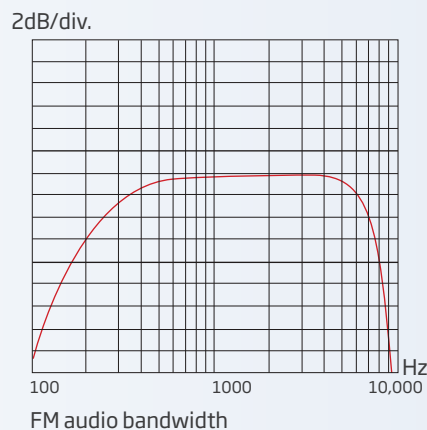

Features that make a difference

ABILITY	FEATURE	WHY THIS IS IMPORTANT
7000 Hz bandwidth	Wide bandwidth	Delivers the broadest signal for natural sound quality and enhanced speech intelligibility.
Auto sleep mode when not in use	Intelligent sleep mode	Allows the child to use Amigo Arc all day long.
Battery and FM LEDs	Visual status indicator	Both caregiver and child can be certain that Amigo Arc is working as it should.
USB charging capabilities	Flexible charging options	Easy, instant charging when necessary.
Remembers previous volume settings	Intelligent volume control	Both caregiver and child can trust that the volume settings will remain comfortable.
Works with Amigo T20/T21/WRP that doubles as wireless synchroniser/programmer	Wireless synchronisation and programming	Synchronisation and programming happens instantly, saving precious time.
Buttons can be programmed to be active or disabled	Tamper resistant controls	Caregivers can be certain the settings are correct.
Low frequency gain emphasis (+6dB)	Superior bass boost	Bass boost in the low frequencies enriches the experience when listening to music.
AUX programmed to IN or OUT	Versatile AUX function	Flexibility to connect either headphones or an external audio source.
Five exchangeable coloured shells	Customisable design	Creates a sense of ownership for children, making them feel more motivated to wear Amigo Arc.

